



Restriction Enzyme

Sda I




Cat.#	Size	Conc.
FG-Sdal	300 units	10 units/μl

Store at -20°C

Supplied with: 10X FastGene® Buffer IV (FG-REB4)
10X FastGene® FastCut Buffer (FG-REBHF)
6X DNA Loading Buffer
Sterile water

Recognition site

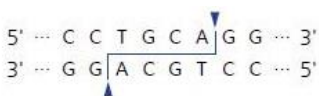

For Research Use Only. Not for use in diagnostic procedures.



Source: *Streptomyces diastaticus* Ng 7-324

Reaction conditions

1X FastGene® IV, 37°C
1X FastGene® FastCut Buffer, 37°C

FastGene® FastCut Buffer

FastGene® restriction enzyme can cut substrate DNA in 5-15 with FastGene® FastCut Buffer.

1X FastGene® Buffer IV

20 mM Tris-acetate (pH 7.9 at 25°C)
50 mM potassium acetate
10 mM magnesium acetate
100 μg/ml BSA

Unit definition

One unit is defined as the amount of enzyme required for complete digestion of 1 μg bacteriophage λ at 37°C for 1 hr in 50 μl reaction mixtures.

Quality control

- Unit definition assay
- Overdigestion assay
- Endonuclease assay
- Extreme pure assay

Standard reaction condition

- Normal protocol

Component	Final Conc.	Volume
Substrate DNA	1 μg	X μl
10X FastGene® Buffer IV	1 X	5 μl
Sda I	10 unit	1 μl
Sterile water		up to 50 μl

→ Incubate at 37°C for 1 hr

- Fast protocol

Component	Final Conc.	Volume
Substrate DNA	1 μg	X μl
10X FastGene® FastCut Buffer	1 X	5 μl
Sda I	10 unit	1 μl
Sterile water		up to 50 μl

→ Incubate at 37°C for 15 min

※ We recommend 5-10 units of enzyme per μg DNA and 10-20 units for genomic DNA in a 1 h digest.

Dilution buffer:

FastGene® Diluent A

Heat Inactivation

Sda I can be inactivated at 80°C for 20 min.

Methylation sensitivity

dam methylation: Not sensitive
dcm methylation: Not sensitive
CpG methylation: Not sensitive

Relative activity in FastGene® Buffers

FastGene® Buffer I:	75%
FastGene® Buffer II:	75%
FastGene® Buffer III:	0%
FastGene® Buffer IV:	100%
FastGene® FastCut Buffer:	100%

Note

It is not affected by *dam*, *dcm*, or mammalian CpG methylation.